


Summary and FAQ for LIS integration



The Sentinel system uses cloud connectivity to provide seamless integration, relying on industry standards AWS and Mirth for robust, high security connection. Our LIS integration will follow a predictable clear process

1. Collaborate to define customized (LIS) Order and Results data
2. Codify and test Order and Results message passing
3. Complete end-to-end testing and verify data syncing

Sentinel Data Flow

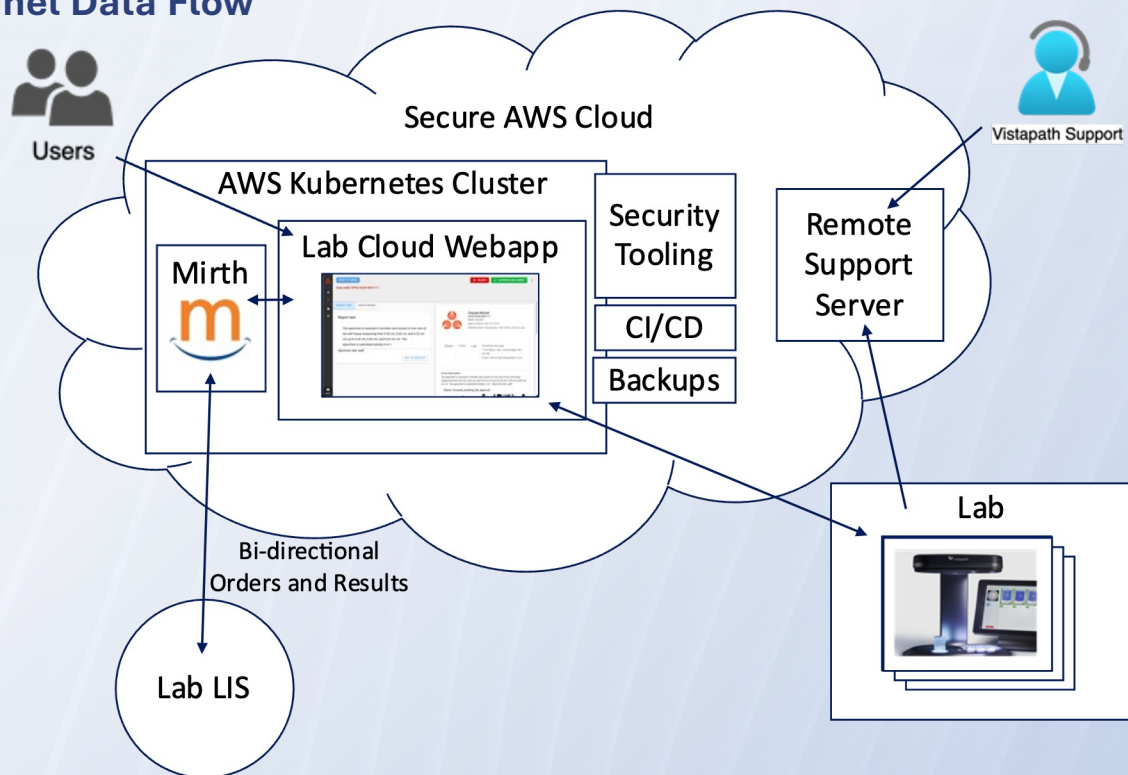

To provide this flexible data flow, Vistapath has designed-in security from the foundation of the Sentinel

- ISO27001 certified
- Data encrypted during transit (TLS) and at rest (AES-256)
- Daily encrypted data backups
- Comprehensive monitoring and alerting with GuardDuty, Cloudtrail, and Falco